


Negative Control

Tube Name: control

Population	#Events	%Parent	control FITC-A Geo Mean	control FITC-A Mean	control FITC-A Median
P1	29,830	35.7	####	353	293
P2	147	0.5	3,519	11,296	2,915
P3	6	0.0	143,265	163,485	161,738

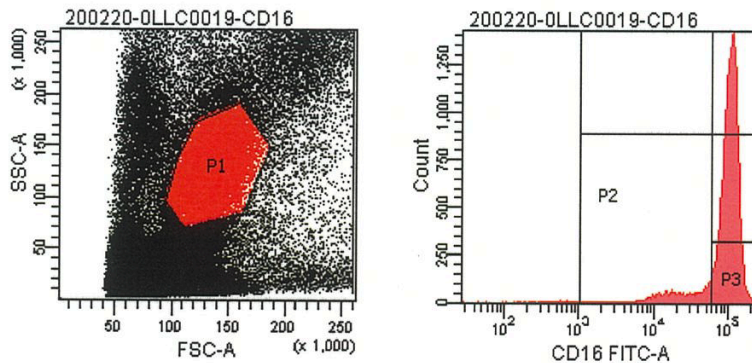

患者サンプル結果

Tube Name: CD16



Population	#Events	%Parent	CD16 FITC-A Geo Mean	CD16 FITC-A Mean	CD16 FITC-A Median
P1	24,798	15.2	####	96,790	102,041
P2	24,615	99.3	85,445	97,506	102,297
P3	21,441	86.5	105,355	108,056	106,809

P3 ゲートのGeo Mean 値で判定、大きく低下する場合はP2ゲート
 (カットオフ値<40000 IGDを強く疑う。<60000 IGDの可能性はある)

報告書2枚目

精度管理

ビーズ: Geo Mean 17000前後に設定

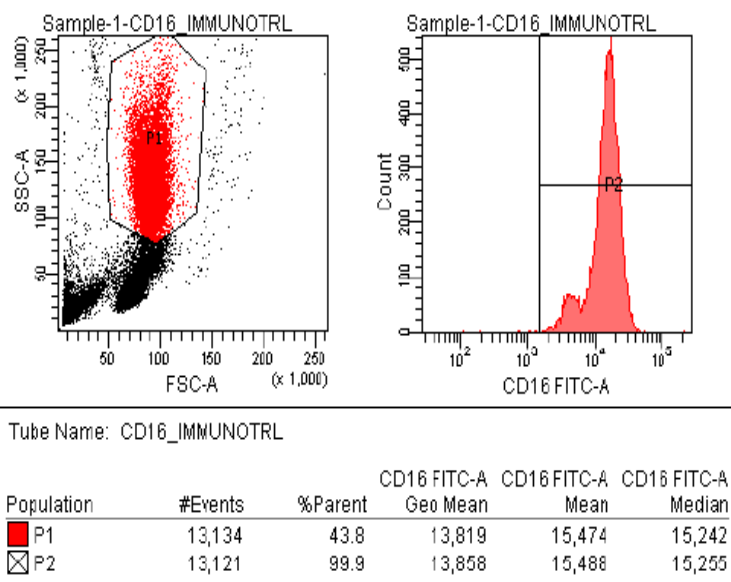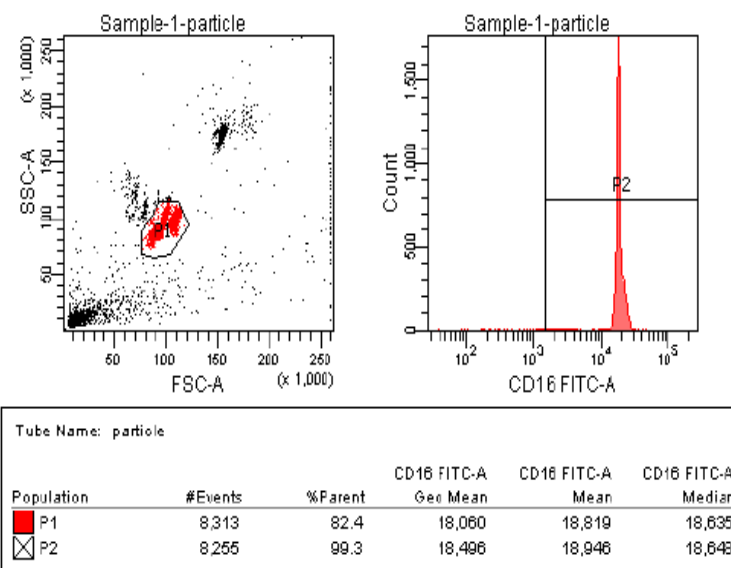

)
:
5

コントロール血球
Geo Mean : 25000前後